

AutoFIX™

6.5 MM TITANIUM CANNULATED

COMPRESSION SCREW

Hexalobe Screw Head design helps to maximize torque while minimizing stripping of the head.

Self Drilling, Self Tapping screws for easy insertion.


Reverse Cutting Flutes for ease of removal and helps minimize screw breakage during removal.

Opposing buttress thread design to optimize the compression.

Fully Cannulated Screw with 1.8 mm K-wire allowing for desired placement of the screw without removing the K-wire.

Strong **Titanium** design to optimize insertion and fixation strength.


Everything needed is in **one easy to use tray**.

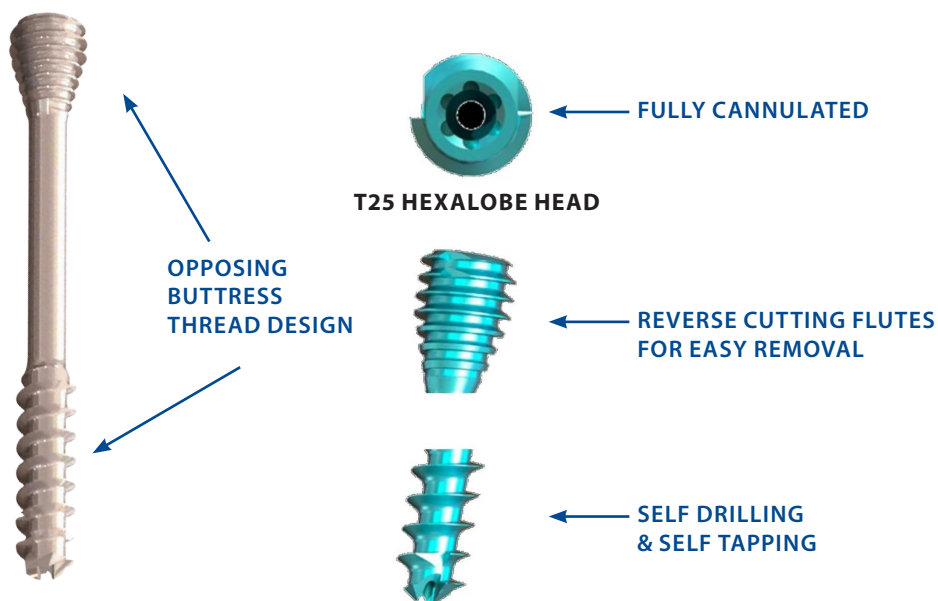

SBi
SMALL BONE INNOVATIONS, INC.



SubTalar Fusion

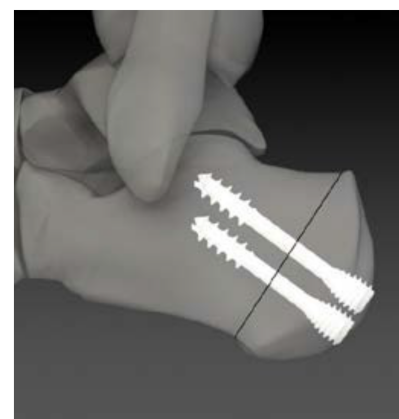


Calcaneal Osteotomy



Ankle Fusion

AUTOFIX 6.5 MM

PART	INSTRUMENT NAME	QTY
241-6545	45mm x 6.5mm	2
241-6550	50mm x 6.5mm	2
241-6555	55mm x 6.5mm	2
241-6560	60mm x 6.5mm	2
241-6565	65mm x 6.5mm	2
241-6570	70mm x 6.5mm	2
241-6575	75mm x 6.5mm	2
241-6580	80mm x 6.5mm	2
241-6585	85mm x 6.5mm	2

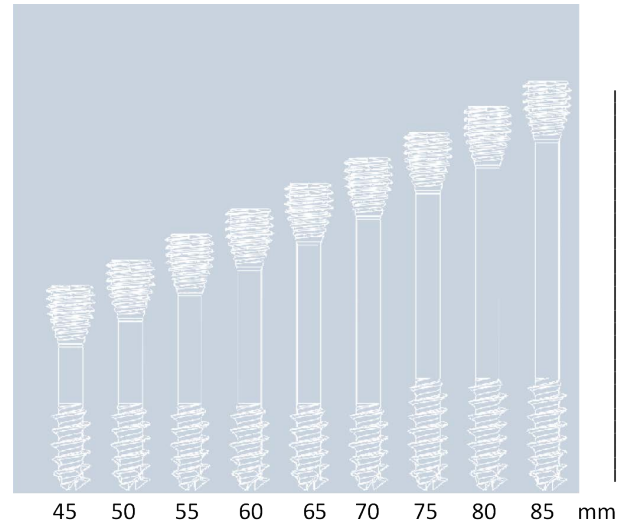

Note: Screw Sizes are available from 45 mm thru 85 mm in 5 mm increments

PART	INSTRUMENT NAME	QTY
603-0026	AutoFIX 6.5 Tray	1
901-0048	Cannulated Screwdriver Shaft T25	2
901-0059	Countersink for AutoFIX 6.5	1
451-0031	K-wire Ø1.8 L150mm	4
909-0022	Depth gauge AutoFIX 6.5	1
903-0054	Cannulated T-handle, AO coupling	1
602-1010	K-wire tube	1
929-0004	Cleaning stylet	1
908-0098	Ø1.8 K-wire Guide	1

Small Bone Innovations, Inc.
 1380 South Pennsylvania Ave.
 Morrisville, PA
 (215) 428-1791/ Fax (215) 428-1795
 SBi Customer Service: (800) 778-8837
 Technical Support: (866) SBI-TIPS

SBI
 SMALL BONE INNOVATIONS, INC.

www.TotalSmallBone.com
www.Star-Ankle.com
www.StarAnkleReplacement.com

Small Bone Innovations International
 ZA Les Bruyères - BP 28
 01960 Péronnas, France
 Tel: +33 (0) 474 21 58 19
 Fax: +33 (0) 474 21 43 12